
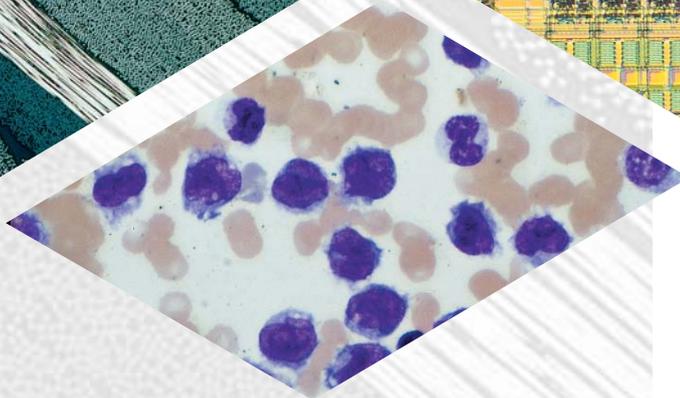
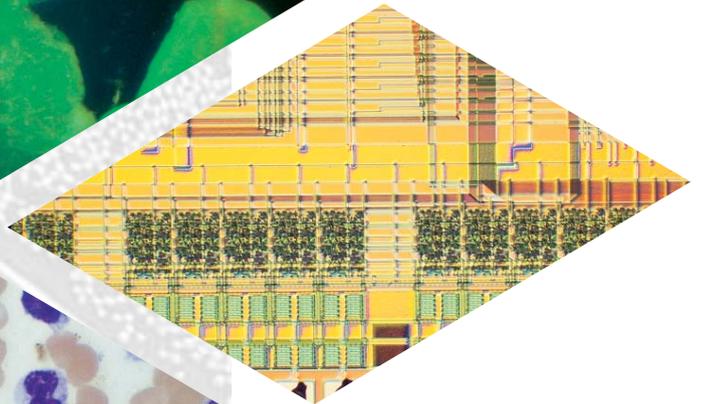
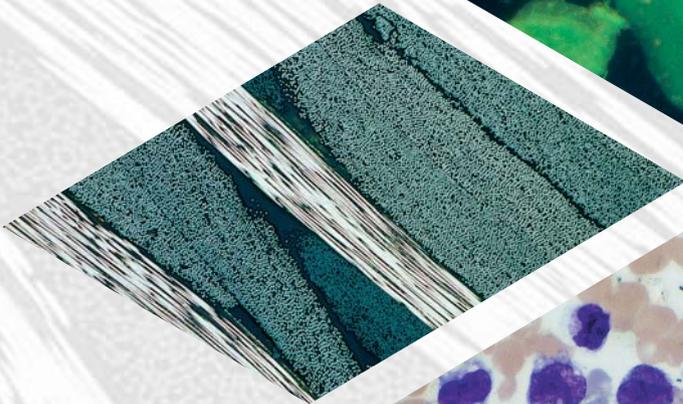


*ProScan™ II*

**PRIOR**  
*Scientific*

**High Performance  
Motorized Stage Systems**

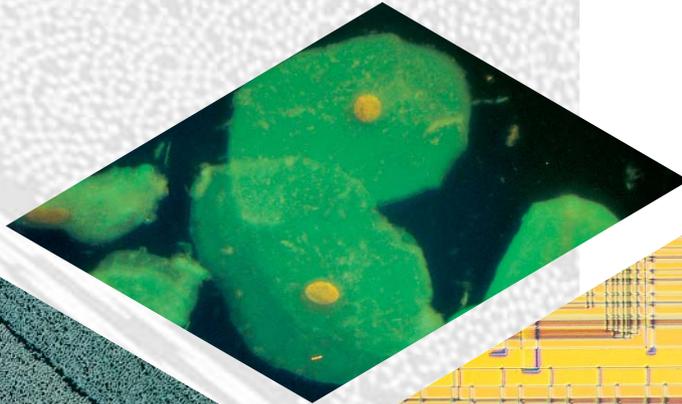

# ProScan II Advanced Microscope Automation



Prior Scientific has been designing and manufacturing precision optical systems, microscopes and related accessories since 1919. This wealth of experience is matched by a commitment to customer service that has earned Prior an enviable reputation for excellent support. These values, plus our understanding of microscopy, provide a unique foundation for the development of an advanced range of motorized stage systems for microscopy and image analysis applications.

Prior proudly introduces the ProScan II system which sets new standards in automated microscopy. Modular by design, a wide range of stages is available for most modern upright and inverted microscopes. Stages may also be adapted to fit other optical inspection systems. The ProScan II controller represents the latest in motor control technology and provides a wide range of advanced features designed for the most demanding applications; not least of which is the capability to control a stage, focussing motor, up to 3 filter wheels and 3

shutters all from a single controller.

Add to this USB and

serial communications,

4 TTL inputs and outputs, encoder feedback for closed loop operation, an optional 4th axis, an advanced autofocus routine and a comprehensive high level software command set including dll's and you can see, the standard is very high indeed. Filter wheels, shutters plus a variety of control accessories such as joysticks and digipots complete the ProScan II product range.

ProScan II is ideally suited to the most demanding imaging applications. A modular design means that systems are easily configured for any combination of stage, focus, filter wheel and shutter options. System performance and reliability are second to none. Furthermore, Prior's flexible, problem solving approach is why the world's leading microscope companies and imaging software manufacturers choose ProScan II.

**HI01**

Stage for upright microscopes providing a travel range of up to 112 x 75mm. Available with a choice of specimen holders, ideally suited to single and four slide applications.

**HI17**

Flat-top stage for inverted microscopes providing a travel range of up to 115 x 77mm. Suitable for slides, multiwell plates, petri dishes, flasks and mounted metallurgical specimens.

**HI05**

Ideal for larger specimens, this stage provides up to 153 x 153mm of travel. Ideal for 150mm semiconductor wafers, photo masks and printed circuit boards.

---

**ProScan II Advanced Controller**


Like the ProScan II range of motorized stages, this advanced controller is designed and manufactured by Prior Scientific. The compact unit enables control of a motorized stage, motorized focus, three filter wheels and three shutters with the speed, accuracy and precision demanded by advanced imaging applications. A fourth axis is also available to enable sample rotation or other custom applications. For triggering camera shutters, relays, linescan cameras or other peripheral devices the unit offers four TTL inputs and outputs. The controller is designed with simplicity in mind. System control using high-level software commands is possible via one of two RS232 ports or the new high speed USB interface. The ProScan II controller has been specifically designed for control by third party software products and the new DLL ensures simple integration. For stand-alone use, without software, the system can be fully programmed and controlled via the Touch Screen Keypad

**Digipot Focus Only Control**

The digipot is ideal for focus only systems. It provides a tactile feel for fine focus adjustments while separate buttons offer immediate control of focus speed and fast movements up and down, for coarse focussing.

**Motorized Focus Control**

ProScan II is ideal for applications where motorized focus control is needed. Step sizes as small as 0.002µm give excellent resolution for precise focus and repeatable positioning in the Z-axis. For large movements when speed is required, the ProScan II focus motor can be driven at speeds of up to 20 revs/s. An optional probe style encoder feedback system provides the highest accuracy and repeatability available.

**Ergonomic Joysticks**

Two or three axis joysticks are available with all stage systems to provide fast, responsive control of the stage. Two programmable Hot Keys are provided which can be used to put a variety of system controls at your fingertips.



**H116**

Provides a travel range of up to 255 x 215mm. The H116 is suitable for large specimens such as 200mm semiconductor wafers and flat panel displays. As with all Prior stages, a range of specimen holders is available.



**H112**

One of Prior's larger stages, the H112 provides up to 302 x 302mm of travel, ideal for 300mm semiconductor wafers.



**HT5050**

One of a range of solid frame stages with X and Y travel of up to 250mm x 150mm respectively, the HT5050 provides 50 x 50mm of travel. Designed for hardness testing machines, the HT stages can withstand loads of up to 100Kg.

with easy to use, menu driven functionality. All system components feature "plug and play" compatibility, allowing the system to automatically detect and configure the attached accessories on start up.

For the most demanding applications ProScan II is compatible with encoders and linear scales to offer the highest repeatability available. These can be fitted to most stages and motorized focus units, after which their output provides a closed-loop capability. This enables use of the system for direct measurement applications where the Touch Screen Keypad becomes a digital read out display.



Other options include autofocus, which uses a unique software algorithm for fast and accurate focussing of specimens. This option is available in both PAL and NTSC formats with six autofocus settings to perfectly match the algorithm with the specific objective lens in use. For specimens that are not uniform the area that the system uses for focus adjustment can be selected from pre-defined quadrants or a user specified focus zone. To ensure that the ProScan II controller is as "future proof" as it possibly can be, the system uses Flash Technology to facilitate software and firmware upgrades.



**Touch Screen Programmable Keypad**

The Programmable Keypad features touch screen technology and a comprehensive yet intuitive software package which provides complete programmability and control of the whole ProScan II system. A range of menus allow for programming of raster, snake and user-defined patterns which can be input, stored and recalled. Information on the current stage position is always available and points of interest can be saved for later review. Patterns and saved points can be downloaded to a PC for further analysis.

**High Speed Filter Wheels and Shutters**

The Filter Wheel system delivers smooth, high speed operation and changes filters in as little as 55ms. Two wheels are available to accept ten 25mm or eight 32mm diameter filters. The filter wheel can be fitted to the excitation and emission ports of your microscope and up to 3 filter wheels and 3 shutters can be operated from 1 controller under PC control or an optional filter wheel keypad. The unique design of this filter wheel allows filters to be changed with ease. A stand to support the illuminator and reduce vibration effects is standard.

An optional sliding filter holder can be fitted for neutral density or infrared filters.



# Precision Motorized Stages

## Features

### Adjustable Limit Switches

Provide the flexibility to reduce the travel range of the stage to match your application and to avoid damaging collisions with your microscope. The limits also provide a precision reference point.

### Wide Range of Specimen Holders

ProScan II stages are available for a variety of applications involving specimens such as slides, multiwell plates, petri dishes, metallurgical samples and semiconductor wafers. Specimen holders are black anodized to provide excellent wear resistance. Custom designs are always considered.

### Linear Encoders

ProScan II stages can be fitted with optional linear encoders for applications requiring greater precision. Encoder output can be fed to a digital read out for measurement applications. In closed loop mode, the encoder output is fed back to the controller to provide the highest repeatability available.

### Precision Ball Screws

High accuracy ground ball screws provide smooth and maintenance free motion. The pre-loaded recirculating ball screw nuts ensure zero backlash. The whole ball screw assembly is connected to the motor with an anti-backlash coupling.

### Cast Aluminium Plate

Prior stages are precision machined out of special cast aluminium plates which are lightweight and provide excellent dimensional stability.

### Precision Stepper Motors


Quiet and precise stepper motors ensure precise positioning of the stage while the use of micro-stepping provides very smooth motion even at low speeds. Motor characteristics offer high acceleration and speeds up to 30 mm/s with 2mm pitch screws.

### Intelligent Scanning Technology (IST)

This new and unique facility allows the controller to electronically interrogate the stage for a range of key performance characteristics and subsequently modify its control output to get the very best stage performance available. The stage model and serial number are also made available to assist with GLP compliance.

### Crossed Roller Bearings


Provide smooth and accurate linear motion for loads of up to 2Kg on open frame stages and 100Kg on solid stages.

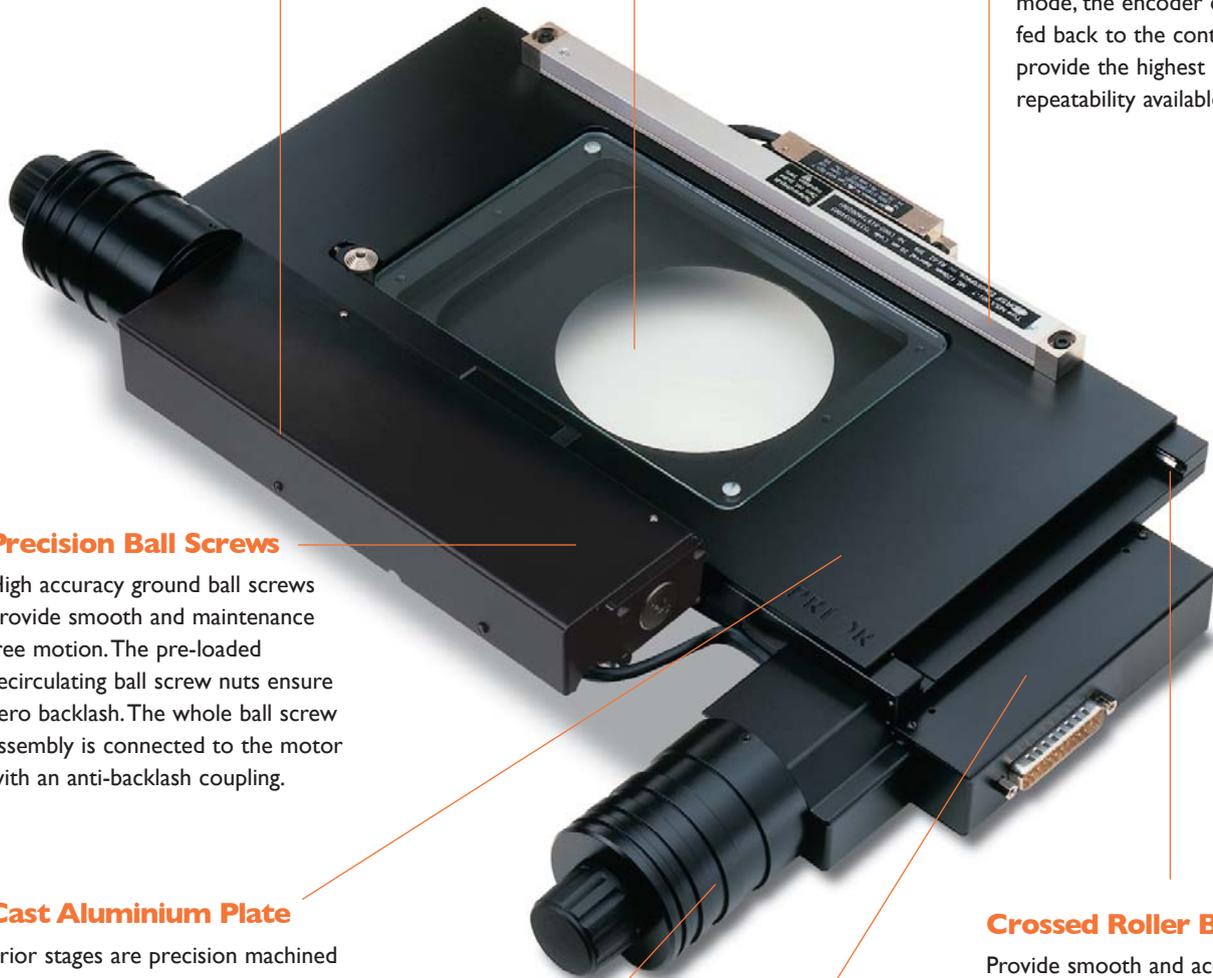

## Specifications

|                       |  |                |  |
|-----------------------|--|----------------|--|
| Power                 | Universal mains input 115/230 VAC, 50-60Hz 200VA                                   | Step Size      | As small as 0.01µm in X and Y and 0.002µm in Z                     |
| Computer Interface    | RS-232C and USB  | Repeatability  | Typically +/- 1µm depending on stage                               |
| RS232 Protocol        | 8 bit word, 1 stop bit, no parity, no handshake, baud rate of 9600, 19200 or 38400 | Linear Slides  | Crossed roller bearings  |
| Controller Dimensions | Width 320mm; Height 90mm; Depth 260mm;   | Drive Screws   | Zero backlash, ground recirculating ball screws, 1, 2 or 5mm pitch |
| Controller Weight     | 2.2Kg (5 lbs)  | Limit Switches | Adjustable in x and y axes (optional in z)                         |
| Stage Speed           | up to 150mm/s (dependant on specifications of stage)                               | Accuracy       | As good as 8µm depending on stage                                  |
|                       |  | Flatness       | 5µm  |

## Specials and OEM Systems

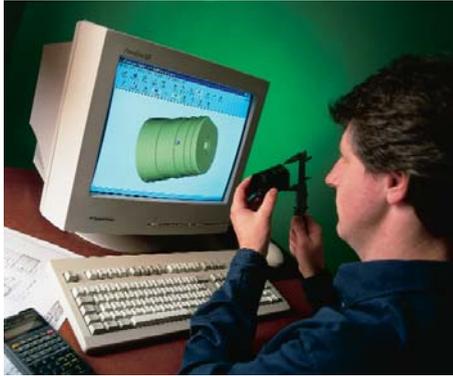

At Prior Scientific we control the design and manufacturing process for all our automated microscopy products. This way, we can be sure of offering the most flexible service. This approach along with our commitment to customer service means that Prior Scientific is uniquely positioned to provide complete systems to match your exact specifications.

The Design Engineering Department employs the latest in Computer Aided 3D modelling along with many years experience in the design and manufacture of scientific instruments. It is here that quality and reliability are designed into our products.

Advanced CNC machines and Computer Aided Manufacture

systems are used to produce high quality components.

In assembly, experienced instrument makers build complete stage and controller assemblies with care and attention to detail.

It is this blend of skills, experience and flexibility that have established Prior as one of the world's leading manufacturers of automated microscopy products. Whether you need a standard product or a custom design, a single unit or OEM quantities Prior Scientific is the right choice!



## ProScan II

### Advanced Microscope Automation



CERTIFICATE NO: FM 61600  
STANDARD: BS EN ISO 9001:2000



**PRIOR SCIENTIFIC INSTRUMENTS LIMITED**  
3-4 FIELDING INDUSTRIAL ESTATE  
WILBRAHAM ROAD FULBOURN CAMBRIDGE CB1 5ET  
TELEPHONE 01223 881711 FAX 01223 881710

**PRIOR SCIENTIFIC INC.**  
80 RESERVOIR PARK DRIVE  
ROCKLAND MA 02370-1062  
TELEPHONE 781-878-8442 FAX 781-878-8736

VISIT PRIOR ON THE WEB AT [www.prior.com](http://www.prior.com)  
Specifications subject to change without notice.